

Kromasil Eternity

Designed for long life


– "The perfect tool for your scale-up activities"

Analytical scale

- 2.5 μm particle size
- UHPLC and HPLC (> 200 000 plates/m)
- Easy scale-up to 5 μm

Preparative scale

- 5 μm particle size
- HPLC
- Easy scale-down to 2.5 μm

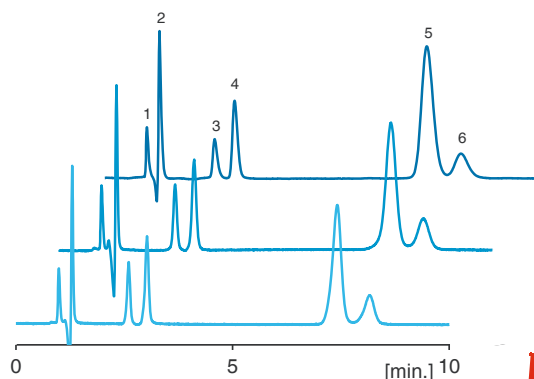

Kromasil Eternity, 5 μm , 21.2 x 50 mm

Kromasil Eternity, 5 μm , 4.6 x 50 mm

Kromasil Eternity, 2.5 μm , 4.6 x 50 mm

Chromatographic conditions

Sample: β -blockers: 1 = Uracil, 2 = atenolol, 3 = pindolol,
4 = metoprolol, 5 = propranolol, 6 = alprenolol
Mobile Phase: acetonitril/water, 50 mM trimethylamine acetate,
pH 11 (40/60)
Flow rate: 0.43 ml/min and 9.0 ml/min
for 4.6 and 21.2 mm i.d. columns, respectively
Temperature: 20 $^{\circ}\text{C}$
Detection: UV @ 230 nm

To order contact
Authorised Distributor:



Hichrom Limited

Tel: 0118 930 3660

Fax: 0118 932 3484

Email: sales@hichrom.co.uk



AkzoNobel